

Viola II

Dynamic Cushion



The Viola II dynamic cushion offers continuous pressure relief and support for seated patients at a Very High Risk of developing pressure ulcers. When combined with an Ultimate Healthcare dynamic mattress system the Viola II offers an affordable preventive package around the clock.

Offering high levels of patient comfort and pressure relief the Viola II can support patients up to 220kg and can be used on a range of chairs.

The compact and lightweight power unit features clear and simple controls that allow the Viola II to be easily set-up for individual patient comfort or clinical requirements, while a stretch, waterproof, vapour permeable cover with an integral covered air pipe and enclosed zip maintain high infection control standards.


Features & Benefits

Alternating Cycle	Moulds to patient shape, providing excellent comfort & pressure re-distribution properties.
Comfort control	Manual pressure adjustment allows the system to be set-up for optimal performance to suit individual patient comfort or clinical requirements.
Audible & visual alarms	Low pressure and power failure alarms enhance patient safety and ensures peace of mind for nursing staff.
Integral foam base & front support (UPRA1717-2 only)	Delivers additional support and security for the patient.
Hanging brackets	Allows the power unit to be attached to a range of chairs.


Specification

Product code:	UPRA1717-2 UPRA2020-2
Risk category:	Up to very high risk
Max user weight:	160kg 25 stone (for UPRA1717-2) 220kg 34 stone (for UPRA2020-2)
Cushion dimensions:	430 x 430 x 110mm (for UPRA1717-2) 500 x 500 x 110mm (for UPRA2020-2)
Power unit dimensions:	145 x 130 x 280mm
Product weight:	2.5kg power unit
Cover material:	Polyester PU
Base material:	Visco elastic foam
Cycle time:	10 minutes
Pressure range:	30 - 70mmHg
Supply voltage:	230V, 50 Hz
Noise level:	60db
Classification:	Class 1 Type BF
Compliance:	EN 60601-1
Warranty:	1 year

Pressure Mapped Patient

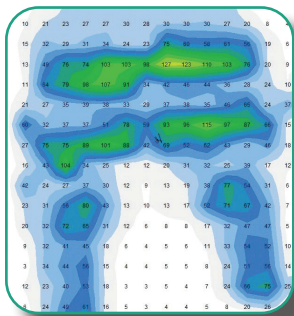

Alternating 1

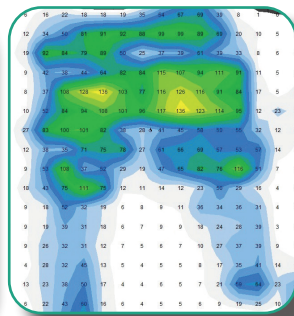

Alternating 2

Age: 61

Gender: Male

Weight: 90kg

Height: 183cm



For more information on this product please contact us on
tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

www.ultimatehealthcare.co.uk