

## Carefree Bariatric

Pressure reducing foam mattress replacement



The **Carefree Bariatric** replaces a standard mattress. It comprises of a foam insert to help provide pressure reduction and features a reinforced core, base and side wall. The mattress provides comfort and protection for the larger patient weighing up to **318kg/ 50 stone**.


The bariatric mattress replacement is ideal for patients considered to be up to a **high risk** of pressure ulcer development and can be used in a **community and long term care** setting.

### Features & Benefits



#### Foam insert

Offering excellent pressure reduction, comfort and support.

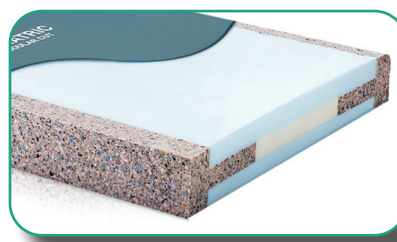

#### Reinforced core, base and side walls

Offering support whilst assisting patients and their caregivers with transfers.



#### Non-turn

Not requiring turning or rotating significant minimising manual handling risks.

#### 2-way stretch, waterproof, vapour permeable PU cover

Improving patient comfort and maintaining skin integrity by reducing shear and friction forces whilst aiding infection control.

For more information on this product please contact us on  
tel: 0333 321 8996 or email: [sales@ultimatehealthcare.co.uk](mailto:sales@ultimatehealthcare.co.uk)

[www.ultimatehealthcare.co.uk](http://www.ultimatehealthcare.co.uk)





## Specification

<b>Product code:</b>	UPRF4878	<b>Dimensions:</b>	122 x 198 x 20cm (48 x 78 x 8")
<b>Max user weight:</b>	318 kg / 50 stone	<b>Fire retardancy:</b>	BS7177 (Crib 5 Composite)
<b>Risk category:</b>	Up to high	<b>Warranty:</b>	2 years
<b>Product weight:</b>	32kg		

**Cleaning:** Can withstand up to 10,000 ppm chlorine (to be used in conjunction with local infection control guidelines) (can be laundered up to 84°C)

## Accessories

### Mattress extension squab infill

<b>USPF480802</b>	for use with	UPRF4878
-------------------	--------------	----------

### Carefree replacement covers

<b>UPRF4878C</b>	122 x 198 x 20cm	for UPRF4878
------------------	------------------	--------------



For more information on this product please contact us on  
tel: 0333 321 8996 or email: [sales@ultimatehealthcare.co.uk](mailto:sales@ultimatehealthcare.co.uk)

[www.ultimatehealthcare.co.uk](http://www.ultimatehealthcare.co.uk)