

Trinity Plus Overlay

Visco elastic foam mattress overlay



The **Trinity Plus Overlay** is used on top of a standard mattress. It comprises of a castellated visco elastic foam to provide pressure redistribution by conforming to a users weight and temperature.

The mattress overlay is ideal for patients considered to be up to **high risk** of pressure ulcer development and can be used in a **community and long term care** setting.

Available in **single** and **double** sizes to provide flexible care options.

Features & Benefits



Castellated visco elastic foam

Allows for immersion, comfort and postural support whilst re-distributing pressure away from high risk areas.



Non-turn mattress overlay

Not requiring turning or rotating significantly minimising manual handling risks.

2-way stretch, waterproof, vapour permeable PU cover

Improving patient comfort and maintaining skin integrity by reducing shear and friction forces whilst aiding infection control.

High frequency welded seams and a fully concealed zip which reduces the possibility of fluid ingress to aid infection control.

For more information on this product please contact us on tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

www.ultimatehealthcare.co.uk






Specification

Product code:	UPRV3478-4	Dimensions:	86 x 198 x 10cm	(34 x 78 x 4") Single
	UPRV5478-4		137 x 198 x 10cm	(54 x 78 x 4") Double
Max user weight:	135 kg / 21 stone		<i>Additional sizes available on request</i>	
Risk category:	Up to high	Fire retardancy:	BS7177 (Crib 5 Composite)	
Product weight:	3kg - 8.5kg	Warranty:	2 years	

Cleaning: Can withstand up to 10,000 ppm chlorine (to be used in conjunction with local infection control guidelines) (can be laundered up to 84°C)

Pressure Mapped Patient

Age: 44	Weight: 83kg
Gender: Male	Height: 170cm

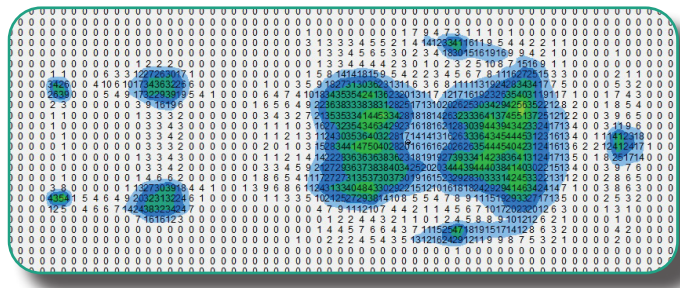

4" Mattress overlay

Accessories

Trinity Plus replacement covers

UPRV3478-4C	86 x 198 x 10cm	for UPRV3478-4
UPRV5478-4C	137 x 198 x 10cm	for UPRV5478-4

Options available

Trinity Plus mattress replacement

UPRV3678	91 x 198 x 15cm	Single (6")
UPRV5478-6	137 x 198 x 15cm	Double (6")



For more information on this product please contact us on tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

www.ultimatehealthcare.co.uk