

# **Trinity Plus Bariatric**

Pressure reducing visco elastic foam mattress replacement



The **Trinity Plus Bariatric** replaces a standard mattress. It comprises of a visco elastic foam insert to provide pressure redistribution by conforming to a users weight and temperature. The mattress also features a reinforced core, base and side wall and provides comfort and protection for the larger patient weighing up to **318kg/50 stone**.


The bariatric mattress replacement is ideal for patients considered to be up to a **very high risk** of pressure ulcer development and can be used in a **community and long term care** setting.



Visco elastic foam insert

Allows for immersion, comfort and postural support whilst re-distributing pressure away from high risk areas.

## Features & Benefits

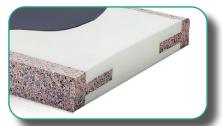

Reinforced side walls, base and core

Offering support whilst assisting patients and their caregivers with transfers.



Non-turn mattress replacement

Not requiring turning or rotating significanting minimising manual handling risks.

#### 2-way stretch, waterproof, vapour permeable PU cover

Improving patient comfort and maintaining skin integrity by reducing shear and friction forces whilst aiding infection control.

High frequency welded seams and a fully concealed zip which reduces the possibility of fluid ingress to aid infection control.



For more information on this product please contact us on tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk



# **Specification**

Product code:	UPRV487801	
Max user weight:	318 kg / 50 stone	
Risk category:	Up to <b>very high</b>	
Product weight:	32kg	

Dimensions:	122 x 198 x 20cm	(48 x 78 x 8")
Fire retardancy:	BS7177 (Crib 5 Composite)	
Warranty:	2 years	

Cleaning: Can withstand up to 10,000 ppm chlorine (to be used in conjunction with local infection control guidelines) (can be laundered up to 84°C)

### **Accessories**

Mattress extension squab infill

**USPF480802** for use with UPRV487801

Carefree replacement covers

**UPRV487801C** 122 x 198 x 20cm for UPRV487801

