

KAKADU

COMMUNE AND SHOWER CHAIR



TABLE OF CONTENT

GENERAL DESCRIPTION.....	2
KAKADU COMMUNE AND SHOWER CHAIR.....	3
FUNCTION AND TECHNICAL DATA.....	4
ACCESSORIES.....	5
KAKADU TILT COMMUNE AND SHOWER CHAIR.....	8

General description

The Mobilex commode and shower chairs are made of high-quality materials, offer the user a high level of comfort and are also easy to use and clean for care personnel.

The commode and shower chairs are height adjustable and can be assembled quickly and easily. When developing our shower and commode wheelchairs, we have focused on making all models easily adaptable so that they can be used as standard for most users.

A wide range of accessories makes it possible to meet most requirements and increase the comfort of the individual user. The stable frame is made of high quality white powder coated aluminium. Screws and other steel parts are consistently made of stainless steel, giving the commode and shower chair a long and maintenance-free life. The seat, back and armrests are made of durable black soft PU foam, which is comfortable for the user and has a timeless design. The four lockable 5" wheels are equipped with an effective protection which ensures that no hair or other dirt can accumulate in the wheel. The wheels are made of plastic and stainless steel and have a doublefunctioning brake.

The Mobilex commode and shower chairs facilitate daily cleaning and can be machine washed in accordance with the latest guidelines, e.g. in home care, nursing homes or similar institutions.

More information?

On these pages you will find a small selection of our extensive product range. If you are interested in our full range, we invite you to visit our website www.mobilex.dk, where you can view all products and download our catalog. The QR code will take you to our website.



KAKADU

COMMUNE AND SHOWER CHAIR



- ① Soft PU foam backrest. A velcro backrest is available as accessory.
- ② Elevating armrests with soft polyurethane foam padded supports.
- ③ Seat with a cut for hygiene. Other types of soft seats can be supplied as accessories.
- ④ Side foldable legrests.
- ⑤ Folding footrest adjustable in height and inclination.
- ⑥ 5" (125 mm) lockable castors, suitable for seat height adjustment.


FUNCTION AND TECHNICAL DATA



The 5" (125 mm) wheels are easy to mount. All four wheels are lockable and equipped with mud flaps.



The legrests can be mounted without tools and can be folded away to the side to facilitate entry and exit.

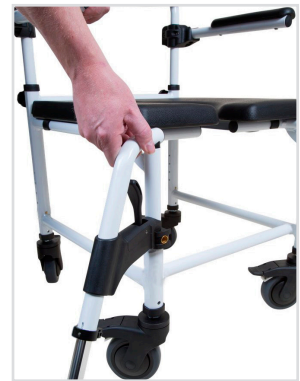

The toilet bucket can be easily removed even when the user is sitting in the chair.



The seat and lid are made of black soft PU foam for comfortable sitting.

Technical data	Small 302015	Standard 302018	Large 302019	XL 302020	Difference with 24"
Seat width	38 cm	45 cm	50 cm	55 cm	---
Seat depth	38 cm	43 cm	43 cm	43 cm	---
Seat height adjust.	50 - 55 cm	50 - 55 cm	50 - 55 cm	50 - 55 cm	---
Seat opening	17 x 33 cm	21,3 x 40,5	21,3 x 40,5	21,3 x 40,5	---
Pass over height **	42,5 - 47,5 cm	43 - 48 cm	42 - 47 cm	40 - 45 cm	---
Pass over width	32 cm	40 cm	45 cm	50 cm	---
Seat/Back angle	90 ° / 12°				---
Distance between armrests	39.5 cm	46,5 cm	51,5 cm	56,6 cm	---
Backrest height	46 cm				---
Total width	49 cm	56 cm	61 cm	66 cm	+ 9,5 cm
Total depth	85 cm	93 cm	93 cm	93 cm	+ 19 cm
Total height adjust.	96-101 cm	96-101 cm	96-101 cm	96-101 cm	---
Max. load	130 Kg	150 Kg	150 Kg	200 kg	---
Weight	11,00 Kg	12,90 Kg	13,50 Kg	14,10 kg	+ 4,20 kg
** +5 cm with the new accessory 30212M					

ACCESSORIES

24" large rear wheels

Item no. 302017

The Mobilex commode and shower chairs can be fitted with large rear wheels so that the user can move around on their own. The 24" rear wheels are supplied complete with parking brakes, pushrims and solid tyres.

The pushrim has an anatomical inner part to provide a good grip even with wet hands. The 24" wheel kit also contains the brakes and all the parts necessary for assembly.

Please note that the 24" wheels increase the total width by 9,5 cm after assembly.



Anti-tip

Item no. 302125

When using the 24" rear wheels, we recommend to use an anti-tip - especially on the "Tilt" model.



Safety belt

Item no. 302122

The safety belt has no metal parts that can rust and is therefore particularly suitable for humid environments. The belt can also be used as an additional support belt for the lumbar region.



Adjustable velcro back

Item no. 302164

To provide the user with optimal back support, a velcro backrest can be adopted which can be optimally adapted to the user's back by means of three adjustable straps. The perforated surface does not absorb water, allowing the skin to breathe.



ACCESSORIES

PU seat and lid

Item no. 302103 (lid) / 302107 (seat)

The seat is made of soft PU foam and is equipped with a hygiene cut-out as standard. The edges of the cut-out are inclined so that the lid sits firmly in it. The seat is fixed with clips and can be easily removed without tools.

The cut-out measures 27 x 19 cm.




PU seat with oval cut-out

Item no. 302109

A soft PU seat with oval cut-out and closed front is optionally available. The seat is attached with clips and can be easily removed without tools.

The cut-out measures 27 x 19 cm.

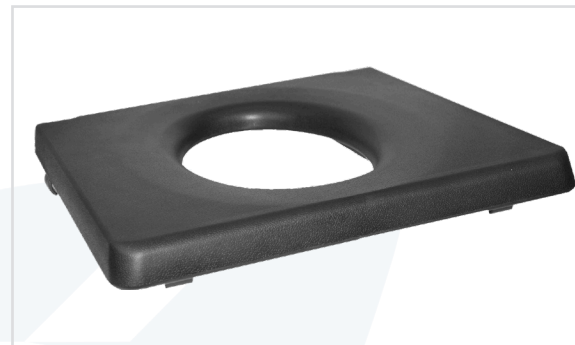

Seat plate lid

Item no. 302110

For both the standard seat and the oval cut-out seat, a closed soft PU seat plate is available which covers the entire seat surface. The seat plate is shaped so that it fits into the seat cut-out and therefore cannot slip, even during movements.



Splash guard

Item no. 30214U

In case the urine flow does not hit the bucket, a splash guard can be attached to the front edge of the hygiene cut-out without tools. The splash guard also acts as a spacer between the legs. The product is made of soft polyurethane and is easy to clean.



Complete head support

Item no. D2101

A headrest is available for all Mobilex commode and shower chairs and can be easily attached to the frame. It is adjustable both in height and depth. The headrest itself is made of soft, comfortable polyurethane foam and can be easily adapted to the shape of the user's head.




ACCESSORIES

Seat upholstery

Item no. 30210U (with open front) / 30211U (closed front)

The commode and shower chairs are equipped with a PU foam seat, which offers good comfort for most people. An additional soft foam cushion covered with Decutex is available for a better pressure relief. All sewings are fully taped so it is easy to clean and 100% waterproof. The cushions are available for both the seat with open front and the seat with closed front.

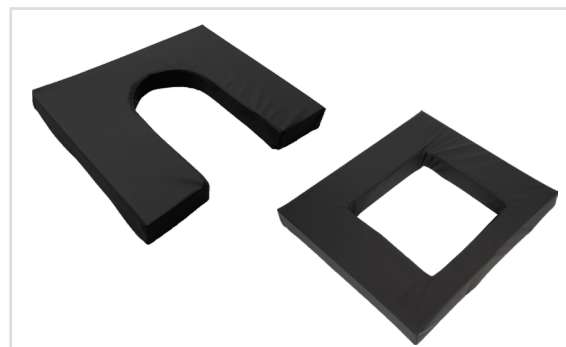

Safety bar

Item no. 302127

In order to offer the user additional safety, a safety bar can be attached to the commode and shower chair. The adapters for the safety bar are fitted to the armrests, after which the bar can be easily removed and fitted.



Amputation support

Item no. 302128 (left) / 302129 (right)

The amputation support supports the leg stump in the correct position and can be fitted in place of the normal footrest without tools. The support is equipped with two moulded polyurethane shells that can be angled independently for optimal support.



Side support pads

Item no. 30212P (small) / 30212Q (large) - left side
30212R (small) / 30212S (large) - right side

The side support pads made of moulded polyurethane are mounted either on the rear tube or on the seat and are infinitely adjustable to the user's needs.

The pads are available in two sizes:
10 x 14 cm
14 x 20 cm



5" wheel with extended tube

Item no. 30212M

There are 5" wheels with extended tube available for the Mobilex commode and shower chairs. This allows the seat height and the passover height of the shower and commode chair to be increased by 5 cm.



KAKADU TILT

The Kakadu Tilt commode and shower chair can be individually adjusted from -5° to +25° in the seat tilt. The headrest with ear recess can be adjusted in height, depth and angle and is included in the delivery.

The lower leg length and the angle of the foot plates are adjustable on the foot rests. The complete black upholstery, including the headrest, is made of soft PUR.

The Kakadu Tilt comes standard with four lockable 5" wheels, toilet bucket, lid and black seat cover.



24" large rear wheels

Item no. 302023

Kakadu Tilt can be fitted with large rear wheels so that the user can move around on their own.

Technical data	Tilt 302024	Tilt L 302025	Diff. with 24"
Seat width	45 cm	50 cm	---
Seat depth	43 cm	43 cm	---
Seat height adjust.	55 - 60 cm	55 - 60 cm	---
Seat opening	21,3 x 40,5	21,3 x 40,5	---
Pass over height **	40 - 45 cm	40 - 45 cm	---
Pass over width	35 cm	40 cm	---
Seat/Back angle	90 ° / 12°		
Distance between armrests	46,5 cm	51,5 cm	
Backrest height	46 cm		
Total width	56 cm	61 cm	+ 14 cm
Total depth	110-120 cm	110-120 cm	+ 19 cm
Total height adjust.	97-141 cm	97-141 cm	---
Max. load	150 Kg	150 Kg	---
Weight	15,10 kg	15,70 kg	+ 4,20 kg
Tilt angle	-5° bis +25°	-5° bis +25°	---

** +5 cm with the new accessory 30212M